

Falcon® Pipet Controller

FALCON®

A Corning Brand

The light, ergonomic, motorized pipet controller is designed for use with glass and plastic serological pipets from 0.5 to 100 mL volume ranges. Conveniently positioned switches change operating modes and speeds quickly to handle different liquid volumes and viscosities. Aspirating and dispensing speed is controlled through the finger triggers. When fully charged, the batteries enable up to 8 hours of continuous use. The large LCD display clearly indicates battery status, pipetting mode and speed.

The Falcon Pipet Controller is supplied with a two-position charging stand, a universal battery charger, and three batteries. Two 0.2 µm and two 0.45 µm hydrophobic polytetrafluoroethylene (PTFE) filters are supplied with the pipet controller. Additional filters are available as standard accessories. The filters, pipet holder and nosepiece are autoclavable.

Features

- ▶ Ergonomic design for exceptional comfort
- ▶ LCD display for continuous visualization of speed, modes and battery status
- ▶ Easy access to the battery compartment
- ▶ Easy maintenance




Falcon® Pipet Controller Ordering Information

Cat. No.	Description
357469	Falcon Pipet Controller, standard version with two 0.2 µm filters and two 0.45 µm filters, 2-position charging stand, universal power supply and set of 3 batteries

Accessories

357472	Hydrophobic PTFE filter (0.2 µm), set of 5
357473	Hydrophobic PTFE filter (0.45 µm), set of 5
357474	Silicone pipet holder
357486	Battery (set of 3)

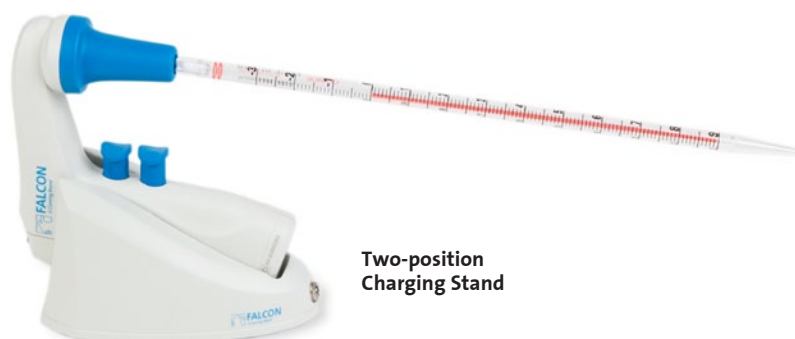

Two-position
Charging Stand

For additional product or technical information, and a complete listing of our International Offices and Distributors, please visit www.corning.com/lifesciences or call 800.492.1110. Outside the United States, please call +1.978.442.2200 or contact your local Corning sales office listed below.

Corning Incorporated Life Sciences

836 North St.
Building 300, Suite 3401
Tewksbury, MA 01876
t 800.492.1110
t 978.442.2200
f 978.442.2476

www.corning.com/lifesciences

Worldwide Support Offices

ASIA/PACIFIC

Australia/New Zealand
t 0402-794-347

China
t 86 21 2215 2888
f 86 21 6215 2988

India
t 91 124 4604000
f 91 124 4604099

Japan
t 81 3-3586 1996
f 81 3-3586 1291

Korea
t 82 2-796-9500
f 82 2-796-9300

Singapore
t 65 6733-6511
f 65 6861-2913

Taiwan
t 886 2-2716-0338
f 886 2-2516-7500

EUROPE

France
t 0800 916 882
f 0800 918 636

Germany
t 0800 101 1153
f 0800 101 2427

The Netherlands
t 31 20 655 79 28
f 31 20 659 76 73

United Kingdom
t 0800 376 8660
f 0800 279 1117

All Other European Countries

t 31 (0) 20 659 60 51
f 31 (0) 20 659 76 73

LATIN AMERICA

Brasil
t (55-11) 3089-7419
f (55-11) 3167-0700

Mexico
t (52-81) 8158-8400
f (52-81) 8313-8589

CORNING | **FALCON®** **cellgro®** **AXYGEN®** **PYREX®** **GOSSSELIN™**

For a listing of trademarks, visit us at www.corning.com/lifesciences/trademarks.
All other trademarks are the property of their respective owners.